


# ROSE K2 SOFT™

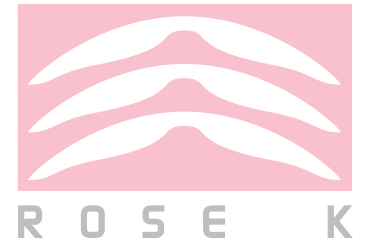

The comfort of a soft lens,  
the simplicity of ROSE K.

## Applications

- **ROSE K2 Soft** is a daily wear soft lens for irregular corneas. ROSE K2 Soft is a 3 month replacement lens manufactured from a silicone hydrogel material.
- **Primary indications:** Intolerance to GP lenses, new contact lens wearers with irregular corneas, early to moderate irregular corneas, if acuity with conventional soft lenses is unsatisfactory, if the environment is unsuitable for GP wear, if a GP lens may be unstable, e.g. sport.

## Design Features

- **ROSE K2 Soft** uses the same simple, proven 5-step fitting method common to all ROSE K designs.
- **ROSE K2 Soft** features an aspheric back optic zone, front surface toricity and front surface aberration control for optimal visual acuity.
- **ROSE K2 Soft** offers precise edge lift control, prism ballast stabilization and reverse geometry for a trouble-free fit.

## Parameter range

- BC range: 7.40 to 9.00 in 0.20 increments
- Diameter range: 14.30 to 15.30mm in 0.10 increments  
Standard diameter: 14.80 mm
- Power range: -30.00D to +30.00D in 0.25D increments
- Cylinder: -0.25D to -10.00D in 0.25D increments
- Cylinder axis: 0° to 180° in 1° increments
- Center Thickness: 0.25mm to 0.60mm in 0.01 steps  
Standard Centre Thickness: 0.35mm
- Edge lifts: 5 Options - Standard (0), Standard Increased (+1.0), Double Increased (+2.0), Standard Decreased (-1.0) and Double Decreased (-2.0)
- Material: Contamac Definitive 74 (efofilcon A): Water content 74%, Dk 60, handling tint

## Diagnostic Set

- Eight lenses in 0.20mm increments from 7.40mm to 8.80mm
- Diameter: 14.80mm
- Center Thickness: 0.35mm
- Material: Contamac Definitive 74 (efofilcon A)
- Standard lift
- Prism ballast, laser marked at base of prism



Customer Care/Ordering: 800.253.9364  
Consultation Direct: 800.566.8001  
Online: [www.artoptical.com](http://www.artoptical.com)